




UV / VIS Spectrophotometer

Single Beam
Double Beam
NANO







Product Overview

	 page 4	 page 5	 page 6
Model	EMC-11S-V	EMC-11S-UV	EMC-11-UV
Single beam	✓	✓	✓
Double beam			
Spectral bandwidth	4 nm	4 nm	4 nm
Wavelength range	325-1000 nm	200-1000 nm	200-1000 nm
Calibration certificate	✓	✓	✓
Cell holder positions	4	4	4
IQ/OQ/PQ	○	○	○
PC Software (Windows®)			
EMC-λ Lambda Software (page 16)	✓	✓	✓
ANALYST software (page 17)			
Audit Trail software (page 18)	○	○	○
Accessories (page 12)			
Peltier/Sipper plug in			
8-position cell changer			
Cell changer up to 50mm	○	○	○
Cell changer up to 100mm	○	○	○
Calibration filter set	○	○	○
Article No.	80 0000 010	80 0000 005	80 0000 016

✓ = included

○ = optional

 page 7	 page 8	 page 9	 page 10
EMC-18S-UV	EMC-61PC-UV	EMC-61PCS-UV	EMC-NANO 2
✓			✓
2 nm	1.8 nm	0.5/1/2/4/5 nm	4 nm
190-1100 nm	190-1100 nm	190-1100 nm	190-1100 nm
✓	✓	✓	✓
4	4+1	4+1	1
○	○	○	○
✓			✓
○	○	○	○
○	○	○	
○	○	○	
○	○	○	○
○	○	○	○
○	○	○	○
80 0000 027	80 0000 031	80 0000 028	80 0000 018

EMCLAB Spectrophotometers are tested with DAkkS (German body of accreditation) certified UV/VIS Reference Materials NIST traceable:

- Photometric Accuracy
- Wavelength Accuracy
- Spectral Resolution
- Stray Light

All EMCLAB Spectrophotometers are supplied with EMCLAB Works Calibration Certificate



Features:

- Color TFT screen
- Self-check system
- Fast choose wavelength
- Auto Zero and Blank
- Sample compartment for different cell holders
- Incl. PC-Software EMC-λ Lambda (Wavelength Measurement, Spectrum Scan, Quantitative Measurement, Kinetic Measurement, DNA/RNA Measurement)
- Incl. EMCLAB Works Calibration Certificate



Our lightest spectrophotometer is especially popular for use in education due to its easy operation and simple functionality. It covers standard photometric and quantitative applications. With the PC software EMC-λ Lambda, kinetic measurements and wavelength scans can also be performed with the instrument.

Model	EMC-11S-V
Wavelength Range	325-1000 nm
Spectral Bandwidth	4 nm
Optical System	Single Beam, grating 1200 lines/mm
Wavelength Accuracy	±2 nm
Wavelength Repeatability	1 nm
Photometric Accuracy	≤±0.5 % T or ±0.004A@1A
Photometric Range	0-200 % T, -0.3 - 3A 0-9999 conc.
Photometric Repeatability	≤0.2 % T
Photometric Mode	T, A, C, F
Stray Light	0.2 % T
Stability	±0.004A/h@500 nm
Baseline Flatness	±0.001A (200-1000 nm)
Noise	0.003A@500 nm
Detector	Silicone Photodiode
Display	Color TFT screen
Central beam height	15 mm
Standard Cell Holder	4-position cell holder 10x10 mm
Light Source	Tungsten lamp
Output	USB port & Parallel port for printer
Power Requirement	AC 85V~265V 50/60 Hz
Dimensions (LxWxH)	440 x 350 x 200 mm
Weight	8 kg
Article No.	80 0000 010

Features:

- Color TFT screen
- Self-check system
- Fast choose wavelength
- Auto Zero and Blank
- Sample compartment for different cell holders
- Incl. PC-Software EMC-λ Lambda (Wavelength Measurement, Spectrum Scan, Quantitative Measurement, Kinetic Measurement, DNA/RNA Measurement)
- Incl. EMCLAB Works Calibration Certificate



The additionally installed deuterium lamp enables the 11S-UV, in contrast to the 11S-V, to also measure within the UV range of the light spectrum. It is therefore suitable for applications that need to cover a greater variability of samples. Kinetic measurements and wavelength scans can also be performed with the instrument via the EMC-λ Lambda PC software.

Model	EMC-11S-UV
Wavelength Range	200-1000 nm
Spectral Bandwidth	4 nm
Optical System	Single Beam, grating 1200 lines/mm
Wavelength Accuracy	±2 nm
Wavelength Repeatability	1 nm
Photometric Accuracy	≤±0.5 % T or ±0.005A@1A
Photometric Range	0-200 % T, -0.3 - 3A 0-9999 conc.
Photometric Repeatability	≤0.2 % T
Photometric Mode	T, A, C, F
Stray Light	0.2 % T
Stability	±0.004A/h@500 nm
Baseline Flatness	±0.001A (200-1000 nm)
Noise	0.003A@500 nm
Detector	Silicone Photodiode
Display	Color TFT screen
Central beam height	15 mm
Standard Cell Holder	4-position cell holder 10x10 mm
Light Source	Tungsten & Deuterium lamp
Output	USB port & Parallel port for printer
Power Requirement	AC 85V~265V 50/60 Hz
Dimensions (LxWxH)	440 x 350 x 200 mm
Weight	9 kg
Article No.	80 0000 005

Features:

- LCD screen (128*64 dots)
- Self-check system
- Auto Zero and Blank
- Sample compartment for different cell holders
- Save the results
- Up to 200 methods & 100 standard curves can be stored
- Incl. PC-Software EMC-λ Lambda (Wavelength Measurement, Spectrum Scan, Quantitative Measurement, Kinetic Measurement, DNA/RNA Measurement)
- Incl. EMCLAB Works Calibration Certificate



The 11-UV is widely used in colleges and enterprises for general quantitative analysis and experiments. The instrument's internal storage offers space for up to 100 standard curves and 200 methods. Via the PC software EMC-λ Lambda, kinetic measurements and wavelength scans can also be performed with the instrument.

Model	EMC-11-UV
Wavelength Range	200-1000 nm
Spectral Bandwidth	4 nm
Optical System	Single Beam, grating 1200 lines/mm
Wavelength Accuracy	±2 nm
Wavelength Repeatability	1 nm
Photometric Accuracy	≤0.5 % T or ±0.003A@1A
Photometric Repeatability	≤0.2 % T
Photometric Range	0-200 % T, -0.3 - 3A, 0-9999 Conc.
Stray Light	0.2 % T
Stability	±0.002A/h@500 nm
Noise	0.003A@500 nm
Detector	Silicone Photodiode
Display	LCD 128*64 dots
Central beam height	15 mm
Standard Cell Holder	4-position cell holder 10x10 mm
Light Source	Tungsten & Deuterium lamp
Output	USB port & Parallel port (printer)
Power Requirement	AC 110/230V 50/60 Hz
Dimensions (LxWxH)	490 x 376 x 220 mm
Weight	14 kg
Article No.	80 0000 016

Features:

- Color TFT screen
- Self-check system
- Auto Zero and Blank
- Sample compartment for different cell holders
- Saves the results
- Up to 200 methods & 100 standard curves can be stored
- Incl. PC-Software EMC-λ Lambda (Wavelength Measurement, Spectrum Scan, Quantitative Measurement, Kinetic Measurement, DNA/RNA Measurement)
- Incl. EMCLAB Works Calibration Certificate



The 18S-UV combines the easy use of the S-series with more comprehensive features of the larger instruments. Wavelength scans and kinetic measurements can be performed directly with the instrument. The instrument's internal storage provides space for up to 100 standard curves and 200 methods.

Model	EMC-18S-UV
Wavelength Range	190-1100 nm
Spectral Bandwidth	2 nm
Optical System	Single Beam, grating 1200 lines/mm
Wavelength Accuracy	±0.5 nm
Wavelength Repeatability	0.3 nm
Photometric Accuracy	≤±0.5 % T or ±0.005A@1A
Photometric Range	0-200 % T, -0.3 - 3A 0-9999 conc.
Photometric Repeatability	≤0.2 % T
Photometric Mode	T, A, C, F
Stray Light	0.05 % T
Stability	±0.001A/h@500 nm
Baseline Flatness	±0.001A (200-1000 nm)
Noise	0.003A@500 nm
Detector	Silicone Photodiode
Display	TFT color screen
Central beam height	15 mm
Standard Cell Holder	4-position cell holder 10x10 mm
Light Source	Tungsten & Deuterium lamp
Output	USB port & Parallel port for printer
Power Requirement	AC 85 V~265 V 50/60 Hz
Dimensions (LxWxH)	490 x 380 x 220 mm
Weight	11 kg
Article No.	80 0000 027

Features:

- Double Beam
- LCD screen (320*240 dots)
- Self-check system
- Auto Zero and Blank
- Sample compartment for different cell holders
- Incl. USB memory stick for methods and results
- Incl. PC software ANALYST (Wavelength Measurement, Wavelength Scan, Quantitative Measurement, Kinetic Measurement, DNA/Protein Measurement)
- Incl. EMCLAB Works Calibration Certificate



With this device, samples and reference are measured simultaneously. This ensures fast and precise measurement results. Variables and values can be entered accurately via the keypad and graphs can be shown directly on the large display. Methods and measurement results can be saved on the instrument's internal storage or an external USB stick.

This instrument features a port for the optionally available Peltier/Sipper System (p. 13).



Model	EMC-61PC-UV
Wavelength Range	190-1100 nm
Spectral Bandwidth	1.8 nm
Optical System	Double Beam, grating 1200 lines/mm
Wavelength Accuracy	±0.3 nm
Wavelength Repeatability	0.2 nm
Photometric Accuracy	≤±0.5 % T or ±0.005A@1A
Photometric Range	0-200 % T, -0.3 - 3A, 0-9999 conc.
Photometric Repeatability	≤0.001A (0-0.5A), ≤0.002A (0-0.5A) ≤0.15 % T (0-100%)
Scan Speed	Hi, MED, LOW, MAX. 3000 nm/min
Stray Light	0.05 % T
Stability	±0.001A/h@500 nm
Baseline Flatness	±0.001A (200-1000 nm)
Noise	0.0003A@500 nm
Detector	Silicone Photodiode
Display	LCD 320*240 dots
Central beam height	15 mm
Standard Cell Holder	4-position cell holder + 1 reference 10x10 mm
Light Source	Tungsten & Deuterium lamp
Output	USB port A for USB memory stick USB port B for PC connectivity Parallel port for printer
Power Requirement	AC 110/230V 50/60 Hz
Dimensions (LxWxH)	589 x 428 x 240 mm
Weight	22 kg
Article No.	80 0000 031

Features:

- Double Beam
- Variable bandwidths 0.5/1/2/4/5 nm
- LCD screen (320*240)
- Self-check system
- Auto Zero and Blank
- Sample compartment for different cell holders
- Incl. USB memory stick for methods and results
- Incl. PC software ANALYST (Wavelength Measurement, Wavelength Scan, Quantitative Measurement, Kinetic Measurement, DNA/Protein Measurement)
- Incl. EMCLAB Works Calibration Certificate



Similar to the 61PC-UV, the 61PCS-UV is a double beam photometer, making it particularly fast and precise. Due to the variable bandwidth, the 61PCS-UV offers maximum freedom in the selection of methods. It is therefore particularly suitable for use in the pharmaceutical sector.

Model	EMC-61PCS-UV
Wavelength Range	190-1100 nm
Spectral Bandwidth	0.5/1/2/4/5 nm
Optical System	Double Beam, grating 1200 lines/mm
Wavelength Accuracy	±0.3 nm
Wavelength Repeatability	0.2 nm
Photometric Accuracy	≤±0.5 % T or ±0.005A@1A
Photometric Range	0-200 % T, -0.3 - 3A, 0-9999 conc.
Photometric Repeatability	≤0.001A (0-0.5A), ≤0.002A (0-0.5A) ≤0.15 % T (0-100%)
Scan Speed	Hi, MED, LOW, MAX. 3000 nm/min
Stray Light	0.05 % T@220, 340 nm
Stability	±0.001A/h@500 nm
Baseline Flatness	±0.001A (200-1000 nm)
Noise	0.0003A@500 nm
Detector	Silicone Photodiode
Display	LCD 320*240
Central beam height	15 mm
Standard Cell Holder	4-position cell holder + 1 reference 10x10 mm
Light Source	Tungsten & Deuterium lamp
Output	USB port A for USB memory stick USB port B for PC connectivity Parallel port for printer
Power Requirement	AC 110/230V 50/60 Hz
Dimensions (LxWxH)	589 x 428 x 240 mm
Weight	24 kg
Article No.	80 0000 028

Features:

- LCD screen (320*240 dots)
- Self-check system
- Auto Zero and Blank
- Unique Flip cell holder 2 in 1 for NANO volume and standard cells
- 0.2 ~ 2.5 μ l sample volume for DNA, RNA, PCR and protein
- Incl. USB memory stick for methods and results
- Incl. PC-Software EMC- λ Lambda
(Wavelength Measurement, Wavelength Scan, Quantitative Measurement, Kinetic Measurement, DNA/RNA Measurement)
- Incl. EMCLAB Works Calibration Certificate



Our new EMC-NANO 2 for measuring samples with nano volumes extends the modern user-friendliness of the 18S-UV with our exclusive nano flip cell holder. This allows both standard samples and nano volumes to be measured easily and accurately. Wavelength scans and kinetic measurements can be performed directly with the instrument.



The unique Flip cell holder for NANO volume and standard cells:

- Simple to use
- NANO cell holder optical path length 0.5 mm
- Cell holder for standard cells optical path length 10 mm
- Horizontal light path - high precision - no loss of energy




Model	EMC-NANO 2
Wavelength Range	190-1100 nm
Spectral Bandwidth	4 nm
Optical System	Single Beam, grating 1200 lines/mm
Wavelength Accuracy	± 0.5 nm
Wavelength Repeatability	0.3 nm
Photometric Accuracy	≤ 0.5 % T or $\pm 0.004A@1A$
Photometric Repeatability	$\leq 0.001A$ (0-0.5A), $\leq 0.002A$ (0-0.5A) ≤ 0.5 % T (0-100%)
Photometric Range	0-200 % T, -0.3 - 3A, 0-9999 Conc.
Scan Speed	Hi, MED, LOW – MAX: 3000nm/min
Stray Light	0.05 % T@220, 340 nm
Stability	$\pm 0.002A/h@500$ nm
Baseline Flatness	$\pm 0.002A$ (200-1000 nm)
Noise	0.0005A@500 nm
Detector	Silicone Photodiode
Display	LCD 320*240 dots
Central beam height	15 mm
Flip Cell Holder	0.2~2.5 μ l and cell 10 mm path length
Light Source	Tungsten & Deuterium lamp
Output	USB port A for USB memory stick USB port B for PC connectivity Parallel port for printer
Power Requirement	AC 85 V~265 V 50/60 Hz
Dimensions (LxWxH)	490 x 370 x 220 mm
Weight	14 kg
Article No.	80 0000 018

UV/VIS Calibration filter set with EMCLAB Works Calibration Certificate for maintenance and service

The EMC-SET-100 for testing the wavelength accuracy and photometric accuracy enables the user to check their measurement results. The EMC-SET-100 is supplied with EMCLAB Works Calibration Certificate of Secondary Spectrometric Calibration Standards.

UV/VIS reference materials for testing:

- Wavelength accuracy
- Photometric accuracy

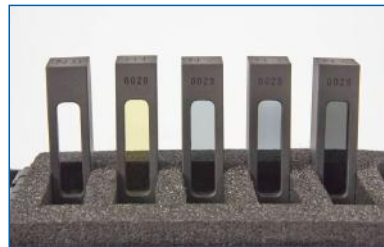
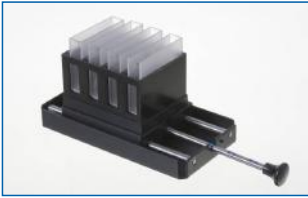



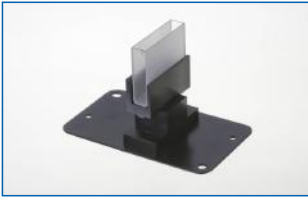



Glass filter set for testing wavelength accuracy (W) and photometric accuracy (A):

Art. No.	Filter	Parameter	Wavelength (nm)
70 0000 045	Holmium Oxide glass filter H1	W*	279; 361; 453; 536; 638 nm
	Neutral density glass filter N2, 0.25 Abs.	A*	440; 465; 546.1; 590; 635 nm
	Neutral density glass filter N3, 0.5 Abs.	A*	440; 465; 546.1; 590; 635 nm
	Neutral density glass filter N4, 1.0 Abs.	A*	440; 465; 546.1; 590; 635 nm
	Empty filter mount N0		

*A= Wavelength for absorbance

*W= Wavelength for wavelength accuracy

Cell Holder	Description	Art. No.
	4-position cell holder for cells 10x10 mm, without slider	70 0000 005
	4-position cell holder for cells up to 50 mm, without slider	70 0000 006
	4-position cell holder for cells up to 100 mm, without slider	70 0000 004
	1-position cell holder for cells 10x10 mm	70 0000 008
	1-position cell holder water-jacketed for cells 10x10 mm	70 0000 009
	4-position cell holder water-jacketed for cells 10x10 mm, without slider (not for EMC-11 series)	70 0000 010
	8-position auto cell changer for cells 10x10 mm (not for EMC-11 series)	70 0000 011
	1-position cell holder for cells up to 100 mm	70 0000 012
	1-position cell holder for solid samples up to 2-5 mm thickness	70 0000 013
	2-position cell holder for solid samples up to 2-5 mm thickness	70 0000 014





Peltier/Sipper System	Description	Art. No.
	Peltier Sipper System EMC-PSA 1, 10°C - 60°C, incl. thermostat-controlled 1-position cell holder and flow through cell 10x10 mm (not for EMC-11 series)	70 0000 086
	Peltier System EMC-PSA 2, 10°C - 60°C, incl. thermostat-controlled 1-position cell holder for cells 10x10 mm (not for EMC-11 series)	70 0000 087
	Sipper System EMC-PSA 3, incl. 1-position cell holder and flow through cell 10x10 mm (not for EMC-11 series)	70 0000 088

Fig. EMC-PSA 1

Lamps	Description	Art. No.
	Halogen (Tungsten) lamp 6V/10W (only for EMC-11S series)	70 0000 078
	Halogen (Tungsten) lamp 12V/20 W (only for EMC-11-UV & EMC-18S-UV)	70 0000 079
	Halogen (Tungsten) lamp 12V/20W (only for EMC-6 series)	70 0000 080
	UV Deuterium lamp type 1	70 0000 081
	UV Deuterium lamp type 2 for EMC-11S-UV	70 0000 082

Cells	Description	Art. No.
	Glass / Quartz glass cells on request	

EMC-11S-V/UV

- A** **Absorbance measurement** Measure the absorbance value of the sample.
- T** **Transmissivity measurement** Measure the transmission of the sample.
- C** **Concentration measurement against a standard sample**
A standard curve through zero is defined by measuring the absorbance value of a standard sample and entering the corresponding concentration. The concentration of the unknown sample is calculated using the standard curve.
- F** **Concentration measurement using the equation coefficients**
A standard curve is defined by entering the coefficients K and B of the standard curve equation $C=K*A+B$, and the standard curve is used to measure the concentration of the unknown sample.
- E** **Energy measurement** For the verification of light source, monochromator and detector.
- System** General system information



Display



Main menu

EMC-11-UV

Basic


Measurements of absorbance, transmission or concentration.
Display and save up to 200 values.

Quantitative

A standard sample with known concentration is used to calculate the concentration of the measured sample.

Utility

Switch lamps on/off, adjust test mode, dark current refresh, wavelength reset, version information

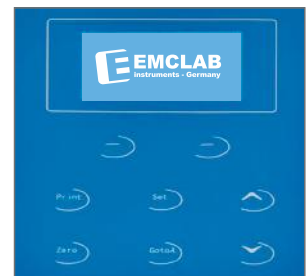

Display



Main menu

EMC-18S-UV and EMC-NANO 2



Photometry Measure the absorbance or transmission value of the sample.



Quantitation The concentration of the unknown sample is calculated using a standard curve.


- a) Standard sample: Measure or enter the absorbance value of a standard sample and enter the corresponding concentration value to define a standard curve.
- b) Equation coefficients: A standard curve is defined by entering the coefficients K and B of the standard curve equation $C=K*A+B$.



Kinetics Measure the photometric value change with time of the sample.



System General system information and settings (e.g. turn lamps on/off, set date and time).

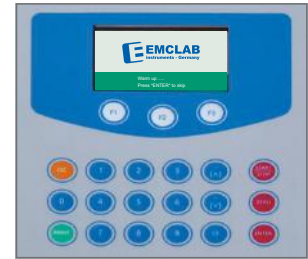

Display



Main menu

EMC-6-series

Main menu: The on-board software includes the following functions: photometry (basic mode), quantitation, wavelength scan, kinetics, DNA/Protein, multiwavelength, system utility.

Multi Wavelength: Measure multiple samples with up to 10 wavelengths.

Spectrum Scan: Choose scan intervals (0.5 to 5 nm), scan speed, as well as the photometric mode to display the spectrum (wavelength/absorbance or wavelength/transmission).

Post processing options include rescaling of axes, curve tracking and determining peaks.

Quantitative Measurement: Define a standard curve using up to 10 standards and choose from 4 curve fitting methods:


1. Linear function
2. Linear function through zero
3. Square function
4. Cubic function

Kinetic Measurement: Absorbance vs time diagrams are displayed in real time. Choose time intervals (0.1 to 60 sec.), delay time as well as the photometric mode to display the curve (absorbance/time or transmission/time).




Post processing options include rescaling of axes, curve tracking and determining the part of the curve required for the rate calculation. The reaction rate is calculated with a linear regression algorithm.

DNA/RNA Measurement: Determine DNA/protein concentration at 260 nm/280 nm or 260 nm/230 nm with optional subtracted absorbance at 320 nm.

Optionally other wavelengths and factors may be entered.

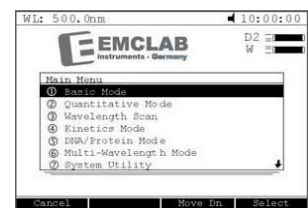

Main menu

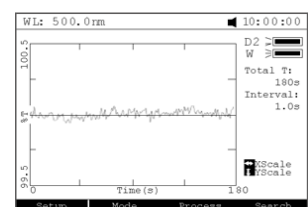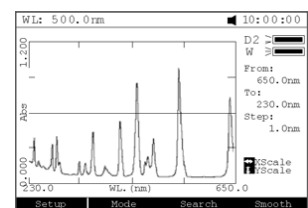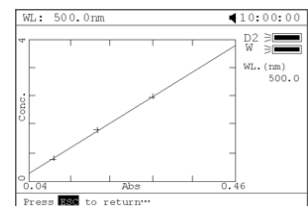

DNA/Protein


No.	Items	Result	Unit
1	A1	0.251	Abs
	A2	0.243	Abs
	Aref	0.095	Abs
	c-DNA	4.524	ug/ul
	C-Pro	110.8	ug/ul
	Ratio	1.059	

PC-Software EMC-λ Lambda for the series EMC-11, EMC-11S, EMC-18S and EMC-NANO 2

The PC-Software is based on Microsoft Windows. EMC-λ Lambda is able to navigate our spectrophotometers via the built-in USB-port. It advances the features of the devices and simplifies the operation by converting your application into easy to set up tasks. The tasks you can choose from include the following:

- Wavelength Measurement**
 With this method it is possible to measure the absorption or transmission of a sample with multiple wavelengths. You may select up to 20 wavelengths and multiple samples per task.
- Spectrum Scan**
 This scan enables you to measure the absorption or transmission of a sample in a specific wavelength range. In accordance with your instrument's specifications, you can choose wavelength range, wavelength intervals, and scan speed.
- Quantitative Measurement**
 This method defines the concentration of a solution in a specific substance. Generate your standard curves by entering parameters or measuring standard samples. To measure your samples, choose any of the saved standard curves from the storage or establish a new one.
- DNA/RNA Measurement**
 With this measurement the concentration for DNA/RNA samples can be measured. You may use common methods or enter parameters to set up a personalized method.
- Kinetic Measurement**
 This method allows you to measure the absorption of a sample at a wavelength over a certain period of time. The results can be used to identify the kinetics of chemical reactions.


Audit Trail PC-Software according to FDA 21 CFR Part 11 optionally available.

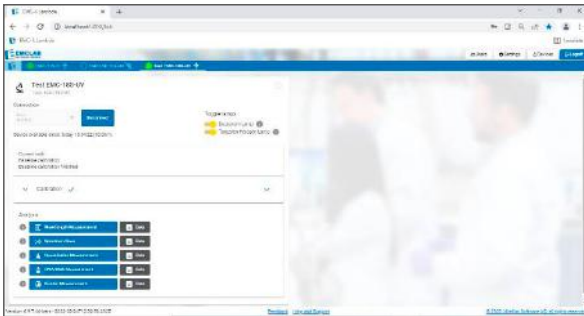


Main Menu

#	wavelength in nm	F2 Abs. @	F3 Abs. @	F4 Abs. @
1	440	0.2639	0.5365	1.0772
2	465	0.2393	0.4877	1.0005
3	546	0.2489	0.5000	1.0068
4	590	0.2868	0.5476	1.0715
5	635	0.2898	0.5348	1.0381

Multi Wavelength Scan: Data

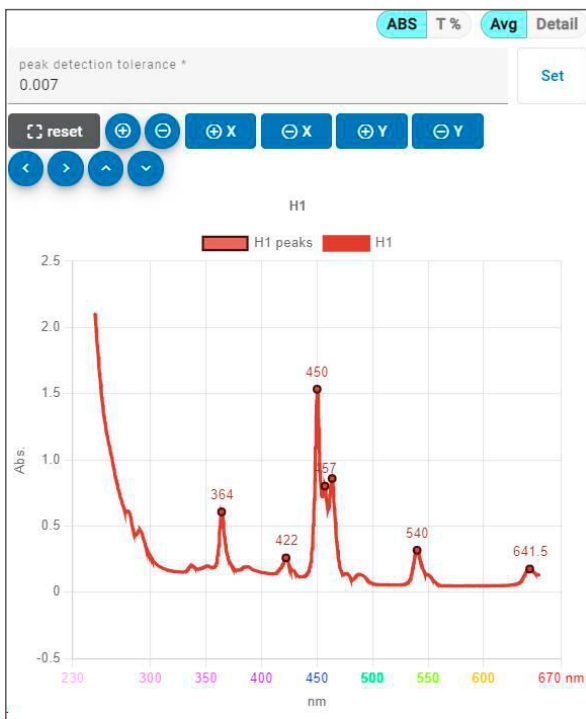

Spectrum Scan: Chart

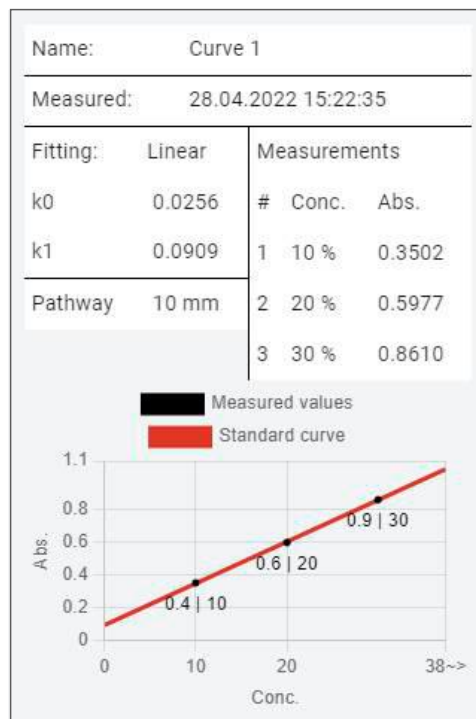

Standard curve

PC software ANALYST for EMC-6 series

The Analyst PC software enhances the functions and data processing and expands the storage capacity. The software includes the following functions:

Wavelength Measurement: Up to 20 different wavelengths can be entered to measure a number of samples.

Spectrum Scan: Choose scan intervals (0.1 to 5 nm), scan speed as well as the photometric mode to display the spectrum (wavelength/absorbance or wavelength/ transmission).

Post processing options include i.a. rescaling of axes, curve smoothing, determining peaks and calculating derivatives.

Quantitative Measurement: Select only one wavelength in the 'method' tab in order to quantitatively calculate the measurement results using a standard curve. Enter or measure up to 20 standards to define a standard curve or enter the function of the standard curve using coefficients. Choose from 4 curve fitting methods:


1. Linear function
2. Linear function through zero
3. Square function
4. Cubic function

Kinetic Measurement: Absorbance vs time diagrams are displayed in real time. Choose time intervals (0.1 to 60 sec.), delay time as well as the photometric mode to display the curve (absorbance/ time or transmission/time).


Post processing options include rescaling of axes, curve tracking and determining the part of the curve required for the rate calculation. The reaction rate is calculated with a linear regression algorithm.

DNA/RNA Measurement: Determination of the DNA/protein concentration of the sample using common methods. Optionally other wavelengths and factors may be entered.


Audit Trail PC software compliant with FDA 21 CFR Part 11 optionally available.

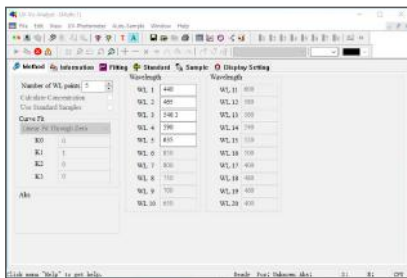

Configuration

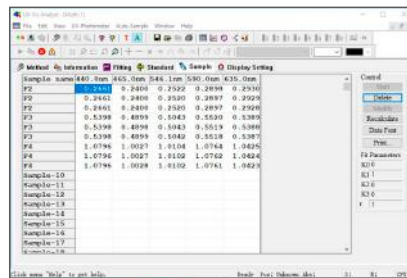

Results

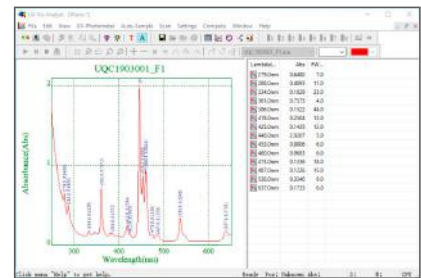

Spectrum Scan

EMCLAB Audit Trail Software (Windows®)

provides secure access of storage results and features compliant with 21 CFR Part 11.

System access requires a user name and password which are assigned by the workgroup manager/administrator. Individual user access levels determine the access to administrative tools which include instrument configuration, analysis applications, user setup, setup and security policies as well as system and application history logs.

For every file version the recorded information is saved (time logged, time created, time last written, deleted, user and manually recorded comments).

EMCLAB Audit Trail is made for

- Monitoring of file modification
- Documentation of changes made by the user
- In specific directories
- Including subdirectories

Article No. 81 0000 010



AUDIT TRAIL

EMCLAB instruments Germany - AUDIT TRAIL

Refresh

Setup

file

- 📁 Data
 - 🚫 Gerät1
 - 🚫 Gerät2
 - 📁 Neuer Ordner
 - 📁 EMC-61PCS-UV
 - 🚫 Neuer Ordner
 - 📁 2018
 - 📁 2017
 - 📁 EMC-NANO-UV
 - 📁 EMC-BR-12

time logged	time created	time last write	length	deleted	comment
📁 file: UQC1707001,UV60_UV600_1.QUA					
07.06.2018 22:45:17	07.06.2018 22:42:23	05.06.2018 15:03:04	1116	<input type="checkbox"/>	
📁 file: UQC1707001_F1_Holmiumoxide_1.sca					
07.06.2018 22:49:29	07.06.2018 22:42:23	07.06.2018 22:49:24	84224	<input checked="" type="checkbox"/>	
07.06.2018 22:45:17	07.06.2018 22:42:23	05.06.2018 15:01:24	84223	<input type="checkbox"/>	
📁 file: UQC1707001_F1-F4.QUA					
07.06.2018 22:50:38	07.06.2018 22:42:23	07.06.2018 22:50:38	1137	<input type="checkbox"/>	Außentemperatur
07.06.2018 22:49:53	07.06.2018 22:42:23	07.06.2018 22:49:53	1137	<input type="checkbox"/>	Zusatz A
07.06.2018 22:45:17	07.06.2018 22:42:23	05.06.2018 15:01:35	1126	<input type="checkbox"/>	Ursprung

Comment

Außentemperatur

K4=0.000000

r=1.000000

[Wavelength]

Num=5

WL1=445.000000

WL2=465.000000

WL3=546.099976

WL4=590.000000

WL5=635.000000

WL6=850.000000

WL7=800.000000

WL8=750.000000

WL9=700.000000

WL10=650.000000

WL11=600.000000

WL12=580.000000

WL13=560.000000

WL14=540.000000

WL15=520.000000

WL16=500.000000

WL17=400.000000


WL18=480.000000


WL19=460.000000

AUDIT TRAIL


Disclaimer



Liquid Handling
 Bottle Top Dispenser
 Titration Station
 Digital Burette
 Digital Pipettes



Stirrer
 Inductive Magnetic Stirrer
 Hotplate Stirrer

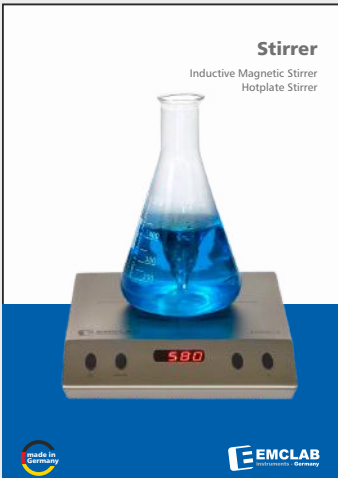

Laboratory Autoclaves
 EMC-ACH 23
 EMC-ACH 29
 EMC-ACH 45







EMCLAB Instruments GmbH
 Kulturstrasse 55
 47055 Duisburg | Germany
 Fon: +49 203 418 59 410
 E-mail: info@emc-lab.de
 Web: www.emc-lab.de

Technical changes reserved